

WHEN YOU PUT PATIENTS' HEALTH FIRST,
REMARKABLE THINGS FOLLOW



State-of-the-art manufacturing facilities.



CHEX
SURGICAL STAPLERS

INFORMED DESIGN.
SCRUPULOUSLY APPLIED.
METICULOUSLY PRODUCED.

OUR RANGE OF PRODUCTS INCLUDE

SINGLE USE CIRCULAR
STAPLER FOR RECTAL
PROLAPSE AND
HAEMORRHOIDS

SINGLE USE
RELOADABLE LINEAR
STAPLER

SINGLE USE CURVED
INTRALUMINAL
CIRCULAR STAPLER

CPH34HV

CS COMPACT™

SINGLE USE
TRANSVERSE CUTTING
LINEAR STAPLER



FOR MORE INFORMATION OR TO VIEW OUR ENTIRE PRODUCT RANGE PLEASE VISIT OUR WEBSITE AT WWW.FRANKENMAN.HK




FRANKENMAN INTERNATIONAL LTD
SUITE B 13/F WING TAT COMMERCIAL BUILDING
121-125 WING LOK STREET SHEUNG WAN HONG KONG
T +852 3106 3035 F +852 3585 0519
E info@frankenman.hk

www.frankenman.hk

Chex and Chex Healthcare are registered
trademarks of Frankenman International Limited



INFORMED DESIGN.
SCRUPULOUSLY APPLIED.
METICULOUSLY PRODUCED.

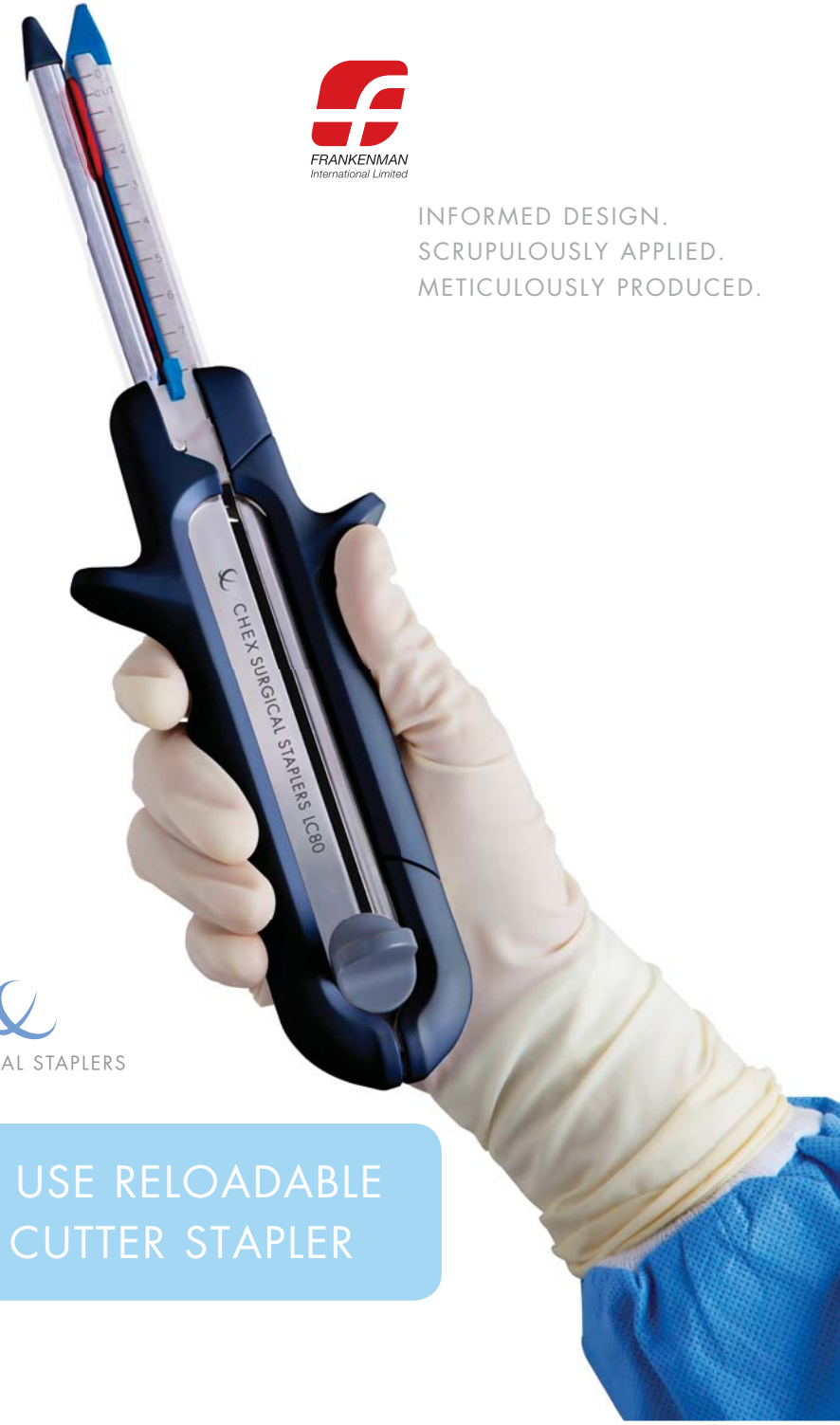

CHEX
SURGICAL STAPLERS

SINGLE USE RELOADABLE
LINEAR CUTTER STAPLER

A RANGE OF PRODUCTS TO MEET THE NEEDS OF MODERN SURGICAL PRACTICE

The Chex™ range of Linear Cutter Staplers focuses on reliability and functionality. Designed to provide accurate and reliable staple formation and tissue resection... time and time again.

At Frankenman, we believe we're not so much in the business of making instruments, but in the business of making surgeries more successful.

That's why we listen to doctors who've put our instruments to use. We read their reports. We invite them to think with us, to brainstorm, to help us develop new and innovative ideas.

That way we know we will be able to deliver the highest quality products to the medical profession both today and way into the future.

THE CHEX™ LC RANGE - KEY FEATURES



SIDE VIEW – ACTUAL SIZE



1
TISSUE RETAINING PIN

2
CLEAR GRADUATIONS AND INDICATION OF 'CUTTING ZONE' ON BOTH SIDES OF STAPLE HOUSING



3
ALL RE-LOADABLE CARTRIDGES HAVE SAFETY LOCK TO PREVENT FIRING OVER EMPTY CARTRIDGE

6
MODERN ERGONOMIC DESIGN DEVELOPED WITH INTERNATIONALLY RENOWNED BRITISH INDUSTRIAL DESIGNERS

7
UNIQUE FIRING 'SAFETY LOCK' TO PREVENT INADVERTENT PRE-FIRING DURING PREPARATION AND POSITIONING



4
UNIQUE 90 DEGREE RELOAD 'SHOULDER' TO ACT AS A POSITIVE TISSUE LIMITER TO AVOID CONCERTINA EFFECT OVER RE-LOAD AND RISK OF PRE-CUT INTERNALLY

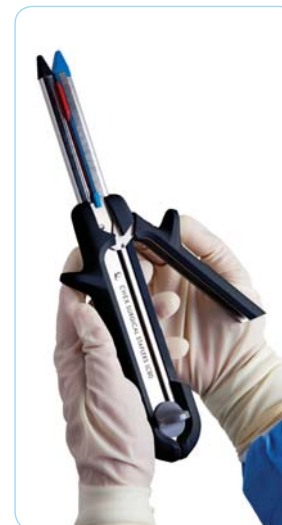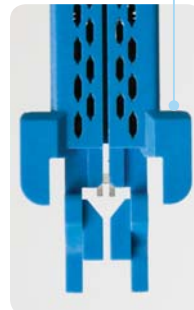

5
CARTRIDGE REMOVAL MADE EASIER. INDENTATION TO EXPOSE CARTRIDGE "LEGS" FACILITATES EASY REMOVAL WHEN CARTRIDGE IS WEDGED BY TISSUE OR FLUID BUILD UP

PRODUCT SPECIFICATIONS

PRODUCT CODE				CLOSED STAPLE SIZE	STAPLING LINE LENGTH	STAPLE QUANTITY
LINEAR CUTTER STAPLER	CARTRIDGE	RELOADS	COLOURED			
LC60G	GREEN	LCR60G	GREEN	2.0MM	61MM	60
LC60B	BLUE	LCR60B	BLUE	1.5MM	61MM	60
LC80G	GREEN	LCR80G	GREEN	2.0MM	81MM	80
LC80B	BLUE	LCR80B	BLUE	1.5MM	81MM	80
LC100G	GREEN	LCR100G	GREEN	2.0MM	100MM	100
LC100B	BLUE	LCR100B	BLUE	1.5MM	100MM	100

KEY FEATURES

- CAM GEAR DESIGN: ENSURES CONSISTENT STAPLE FORMATION HELPING PREVENT BLEEDING.
- EMPTY CARTRIDGE SAFETY LOCK: CUTTING BLADE WILL ONLY ACTIVATE WITH FULL CARTRIDGE.
- RELOADABLE COLOUR CODED CARTRIDGES.